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***BRUGSANVISNING***

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# USER MANUAL

## Inflatable bathtub with electric pump

Product code: EASY BATH

### 1. Characteristics

Portable, inflatable bathtub is made of durable vinyl. Dual air chambers provide a protection against spillage of water and suds. Water drain allows for quick and safe emptying of bathtub after usage.

**The set includes:** bathtub, electric pump, repair kit

### 2. Application

The inflatable bathtub allows the patient's body to be washed in the supine position. It is particularly recommended for indoor use during care for the bedridden and elderly people with difficulties in tilting their head and bending.

### 3. Method of use

**Step 1.** Put the patient inside the bathtub. Inflate the bathtub using the pump and close the water drain valve at the bottom of the bathtub. The bathtub can be filled with warm water with temperature below 50°C.

**Step 2.** After the patient is, open the water drain valve to empty the bathtub through the drain.

**Step 3.** Release the air through the pump valve and transfer the patient. Dry the bathtub.

#### 4. Technical data

|                    |                  |
|--------------------|------------------|
| Bathtub dimensions | 190 x 90 x 22 cm |
| Bathtub weight     | 2,00 kg          |
| Pump weight        | 0,20 kg          |
| Bathtub volume     | 35 liters        |
| Pump power usage   | 150W             |

#### 5. Notes

**Timago International Group does not bear responsibility of improper use of the bathtub, failure to comply safety regulations and misuse.**

#### 6. Maintenance

- Clean all parts of the set after use and allow them to dry.
- When not in use, protect from high temperature and in a humid environment.
- Do not expose the kit to direct contact with sharp objects.
- Adjust water temperature to provide comfort for the patient.
- Do not use detergents, solvents or essential oils.
- Do not expose the electric pump to direct contact with water.
- The pump operation mode should not last longer than 20 minutes with 15 minutes break. **Do not use pump while using the bathtub.**

#### 7. Storage and transport conditions

Device should not be exposed to direct sunlight, water or moisture.

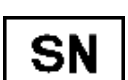
#### 8. Warranty information

All products distributed by our company are covered by a warranty, the terms of which are described in the warranty card available on our

website. Please contact the retailer from whom you purchased the product. Please, keep in mind that for warranty purposes, the proof of purchase (receipt or invoice) should be retained.

## 9. Labels

|                                                                                    |                                                                                      |                                                                                   |                             |
|------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------|-----------------------------|
|   | Reference number                                                                     |  | Manufacturer                |
|   | LOT number                                                                           |  | Manufacture date            |
|   | Serial number                                                                        |  | Medical device              |
|   | Note                                                                                 |  | Please read the instruction |
|  | The manufacturer checked for compliance with essential medical devices requirements. |                                                                                   |                             |

Our complete line of products is available from our distributors.

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